



marClamp[®] IQ

marClamp[®] IQ Instruments and maXium[®]



The bipolar vessel sealing system
SealSafe® IQ for open surgery

marClamp® IQ Instruments for bipolar vessel sealing with SealSafe® IQ

The marClamp® IQ reusable instruments are used in combination with the maxium® HF generator for bipolar vessel sealing and enable effective, lasting and safe sealing of vessels and tissue bundles in all open surgical applications.

marClamp® IQ instruments can be combined with all maxium® units from hardware 06.

“Plug and Play” – automatic recognition of instruments with smart IQ plugs

marClamp® IQ instruments are equipped with a special plug by which the instrument is automatically assigned the required type of current, SealSafe® IQ, and an output power suitable for the instrument.

Intelligent IQ technology makes it possible to use marClamp® IQ instruments both in “Plug and Play” mode and by means of the accustomed program philosophy.

Intelligent vessel sealing with marClamp® IQ and SealSafe® IQ

maxium®		
+ marClamp® IQ		
= SealSafe® IQ		
... the more intelligent sealing system!		

SealSafe® IQ

The bipolar vessel sealing system consisting of the maxium® high-frequency generator with its SealSafe® IQ current and the bipolar marClamp® IQ instruments allows you to seal vessels and tissue bundles effectively and lastingly without any need for previous dissection or detailed exposure of the tissue to be sealed.



Thanks to the precise SealSafe® IQ current specially adapted to this kind of application, only the tissue located between the instrument's jaws is sealed – lateral thermal tissue damage can thus be kept to an absolute minimum. This is explained by the high current-to-voltage ratio (high current but low voltage) that is typical of SealSafe® IQ current, distinguishing it from the conventional bipolar HF technique.

Besides, the SealSafe® IQ current has been optimized to reduce tissue adhesion or instrument incrustation drastically and prevent tissue carbonization as long as the current is used as intended.

Due to the effective combination of high pressure (p) and electric current (I), a permanently sealed zone is reliably created on arteries, veins or tissue bundles. The instrument's ratchet plays a vital role here because it is important for the pressure to be maintained and kept constant throughout the sealing process.

During this process, endogenous structural proteins such as collagen and elastin are transformed so that a permanently sealed zone is created as a result.

With this sealing method there is absolutely no need for exogenous materials such as sutures or clips. The advantage of this method is the reduction in post-operative pain because no changes in tissue structures occur, e.g. constriction with suture material.

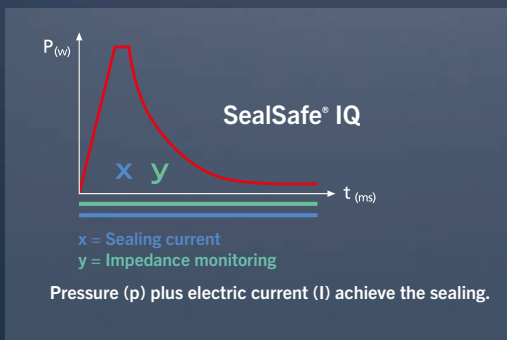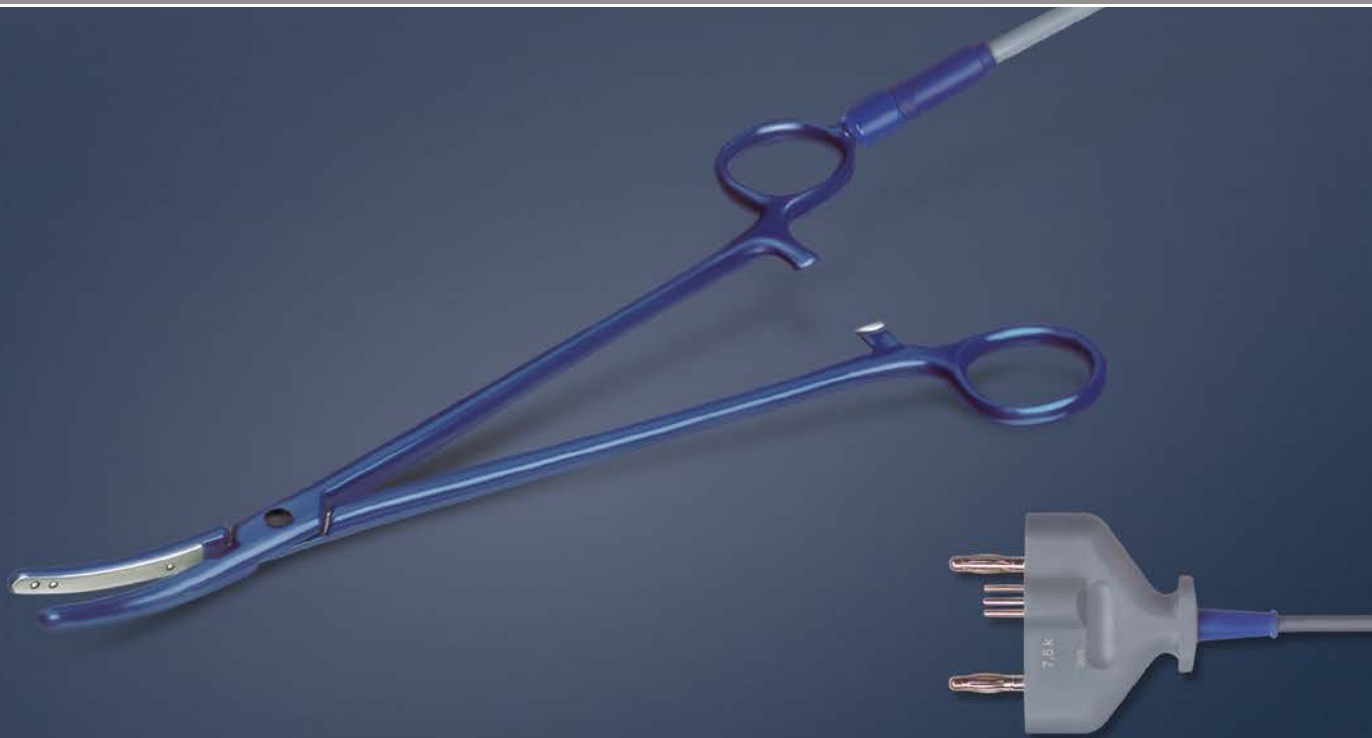


marClamp® IQ applications

- General and visceral surgery
 - Gastrectomy
 - Sigmoidectomy
 - Colectomy
 - Cholecystectomy
 - Mastectomy
 - Splenectomy
- Gynecology
 - Abdominal hysterectomy
 - Vaginal hysterectomy
 - Adhesiolysis
 - Cystectomy
- Urology
 - Prostatectomy
 - Cystectomy
- ENT
 - Thyroidectomy
 - Strumectomy
 - Parotidectomy
 - Neck dissection
- Proctology
 - Hemorrhoidectomy



Advantages of the system

- Sealing process for vessels up to 7 mm with a burst pressure of over 360 mm Hg
- Minimal lateral heat spread of less than 2 mm
- Reliable sealing of vessels or tissue bundles without any previous dissection and detailed exposure
- Precise application due to limited lateral heat spread
- Less encrustation and carbonization thanks to the maxium® generator's special SealSafe® current
- Reliable surgical results without any need to use exogenous materials such as suture or clips
- Time savings compared with suturing; a single application is usually enough → shorter operating times
- Autoclavable at 134°C, acc. to ISO 17665-1
- High efficiency with reusable instruments, no expenses for consumables
- Instruments with ratchet for defined pressure application
- Ceramic pins in one jaw of the instrument ensure a defined distance between the jaws and prevent tissue slippage
- "Plug and Play" function with IQ plug (from hardware 06)
- Compatible with existing KLS Martin standard bipolar connecting cables
- Simple, yet highly reliable design
- No extra costs if the multidisciplinary HF unit maxium® is used / no special equipment required

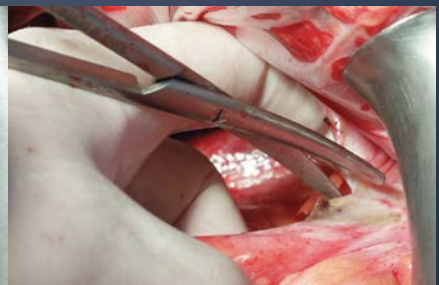

The sealing process is started by grasping the tissue with the bipolar marClamp® IQ instrument and closing the ratchet. This strongly compresses the tissue.

The HF current is now activated via the foot switch until the sealing process stops automatically. Completion is confirmed by an acoustic signal.

Finally, the sealed tissue can be safely transected with a conventional pair of scissors.



The instrument can be reused after autoclaving it at 134°C (273.2°F). Due to the special current flow, no insulation is needed in the joint or box lock. This makes it highly resistant to wear.

One of the instrument's jaws features ceramic pins that ensure a defined distance between the jaws over their entire length, thus preventing a short. Besides, the pins prevent tissue slippage as well.

To apply the necessary pressure to the tissue and maintain it throughout the sealing process, the marClamp® instruments feature a ratchet.

Overview of marClamp® IQ models

Different instrument geometries give each user the right instrument for the task at hand.

marClamp® IQ
maXium®





Please note that the marClamp® IQ models are intended for exclusive use with the maxium® HF generator and its SealSafe® IQ current.



marClamp® IQ

- ✓ Plug & Play – smart bipolar vessel sealing forceps with instrument identification function
- ✓ Automatic parameterization on the maxium® HF generator
- ✓ Fully complies with the proven maxium® program philosophy

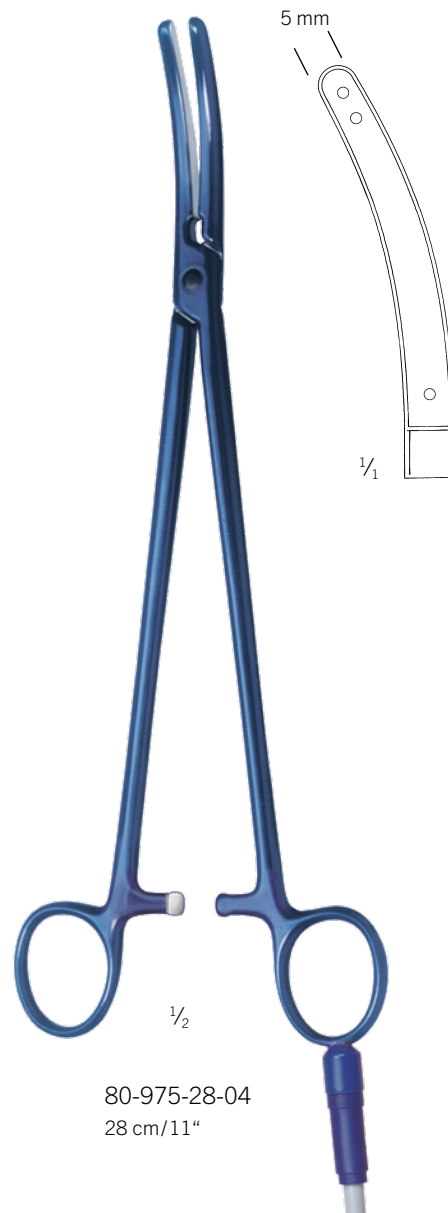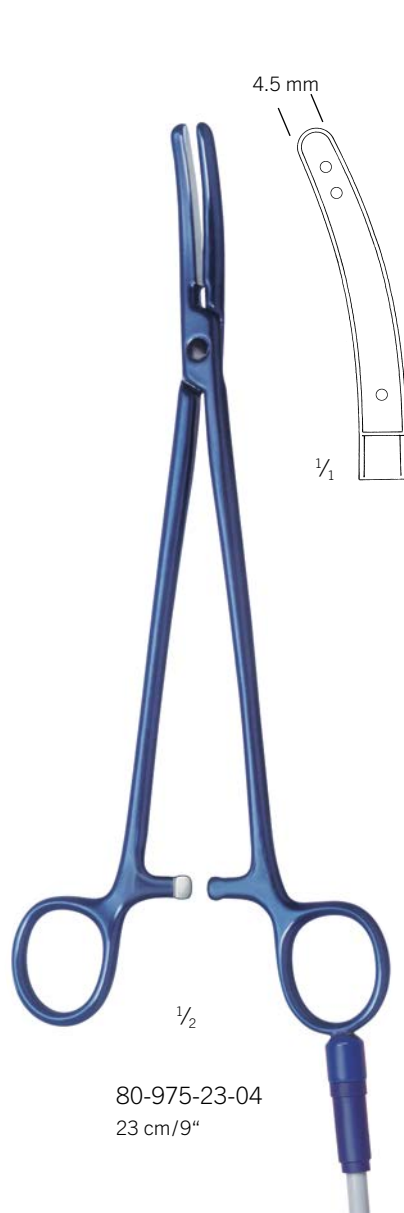



Overview of marClamp® IQ models

Different instrument geometries give each user the right instrument for the task at hand.

marClamp® IQ
maXium®





Please note that the marClamp® IQ models are intended for exclusive use with the maxium® HF generator and its SealSafe® IQ current.

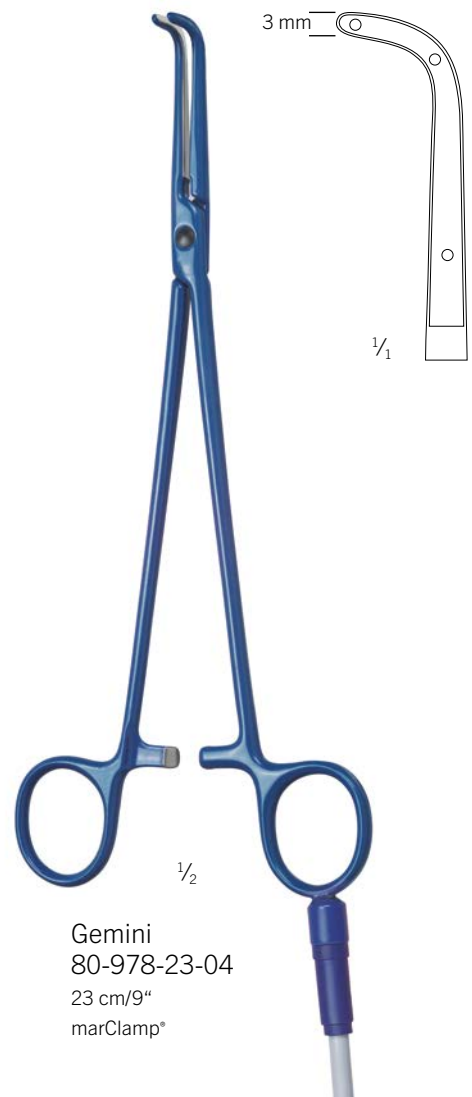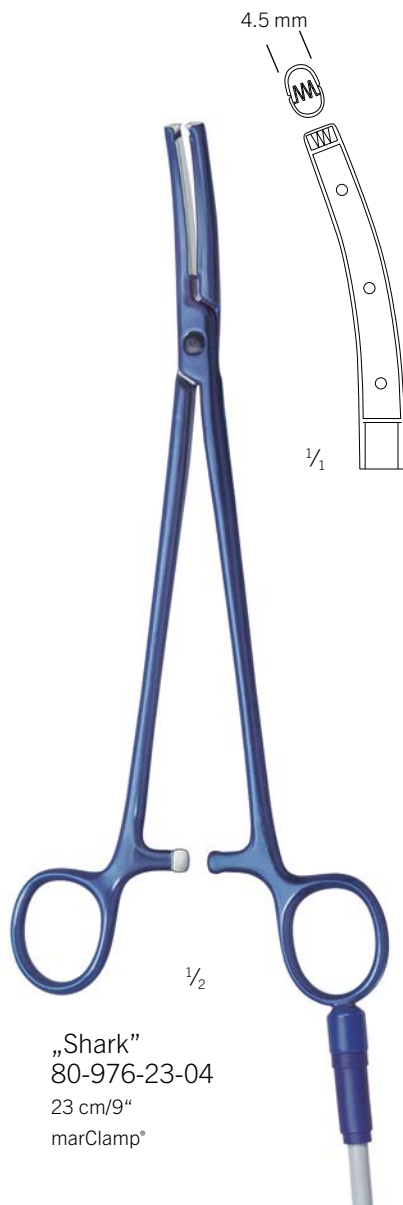

marClamp® IQ

- ✓ Plug & Play – smart bipolar vessel sealing forceps with instrument identification function
- ✓ Automatic parameterization on the maxium® HF generator
- ✓ Fully complies with the proven maxium® program philosophy

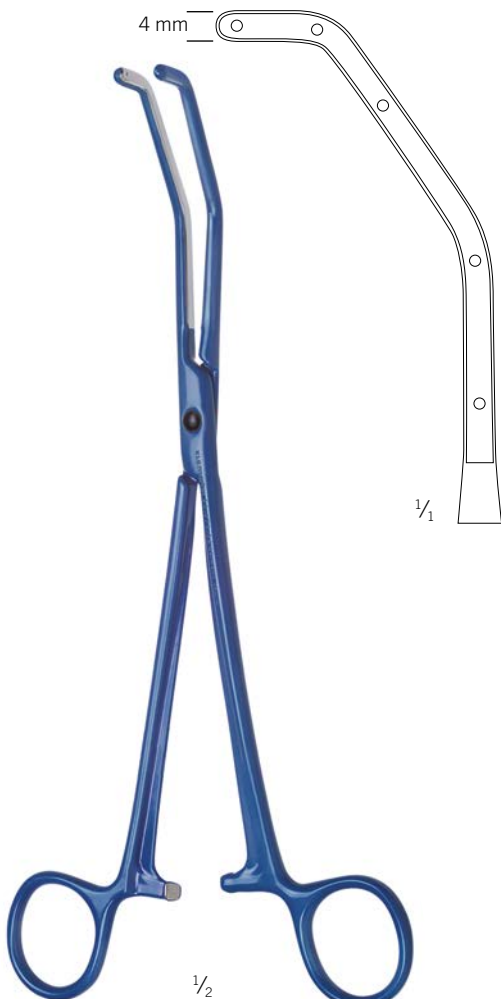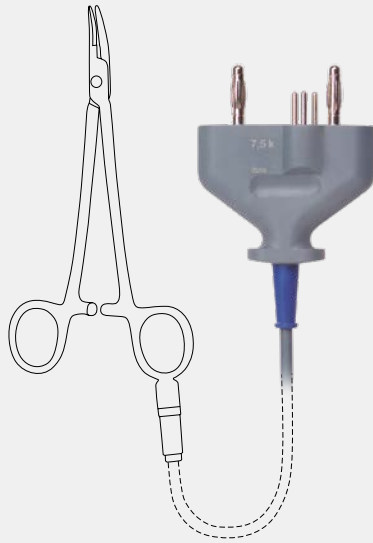

Kalem
80-977-23-04
25 cm/9 7/8"
marClamp® Uro

KLS Martin Group

KLS Martin Australia Pty Ltd.

Sydney · Australia
Tel. +61 2 9439 5316
australia@klsmartin.com

KLS Martin do Brasil Ltda.

São Paulo · Brazil
Tel. +55 11 3554 2299
brazil@klsmartin.com

KLS Martin Medical (Shanghai) International Trading Co., Ltd

Shanghai · China
Tel. +86 21 5820 6251
info@klsmartin.com

KLS Martin India Pvt Ltd.

Chennai · India
Tel. +91 44 66 442 300
india@klsmartin.com

KLS Martin Italia S.r.l.

Milan · Italy
Tel. +39 039 605 67 31
info@klsmartin.com

KLS Martin Japan K.K.

Tokyo · Japan
Tel. +81 3 3814 1431
info@klsmartin.com

KLS Martin SE Asia Sdn. Bhd.

Penang · Malaysia
Tel. +604 261 7060
malaysia@klsmartin.com

KLS Martin de México, S.A. de C.V.

Mexico City · Mexico
Tel. +52 55 7572 0944
mexico@klsmartin.com

KLS Martin Nederland B.V.

Huizen · Netherlands
Tel. +31 35 523 45 38
infonl@klsmartin.com

KLS Martin SE & Co. KG

Moscow · Russia
Tel. +7 499 792 76 19
russia@klsmartin.com

KLS Martin Taiwan Ltd.

Taipei · Taiwan
Tel. +886 2 2325 3169
taiwan@klsmartin.com

KLS Martin SE & Co. KG

Dubai · United Arab Emirates
Tel. +971 4 454 16 55
middleeast@klsmartin.com

KLS Martin UK Ltd.

Reading · United Kingdom
Tel. +44 118 467 1500
info.uk@klsmartin.com

KLS Martin LP

Jacksonville · Florida, USA
Tel. +1 904 641 77 46
usa@klsmartin.com

KLS Martin SE Asia Sdn. Bhd.

Hanoi · Vietnam
Tel. +49 7461 706-0
info@klsmartin.com

Do you know how to get all important information about the KLS Martin Energy Devices?

Please download KLS Martin App Energy Devices! The App is available for Android and iOS.



<https://itunes.apple.com/de/app/klsmartin-energy-devices/id1198171415?l=de&ls=1&mt=8>



<https://play.google.com/store/apps/details?id=com.klsmartin.energydevices>



KLS Martin SE & Co. KG

A company of the KLS Martin Group

KLS Martin Platz 1 · 78532 Tuttlingen · Germany
PO Box 60 · 78501 Tuttlingen · Germany
Tel. +49 7461 706-0 · Fax +49 7461 706-193
info@klsmartin.com · www.klsmartin.com