

smartPen

The new monopolar electrode handpieces
Touch and feel the difference!

smartPen

Premium quality electrode handpieces for your daily work in the OR

HF electrosurgery is indispensable in daily routine OR work. Apart from bipolar electrosurgery, e.g., with bipolar forceps, combined with an active electrode monopolar handpieces are used to cut, dissect, and coagulate.

With the new generation of electrode handpieces KLS Martin offers the surgeon a device designed to fulfill today's ergonomic demands.

- ✓ Ergonomic
- ✓ Premium quality
- ✓ Safe

Overview smartPen electrode handpieces

Item No.	Electrode holder	Activation	Plug
Electrode handpieces with permanent cord			
80-217-13-04	4.0 mm	2 buttons	KLS Martin
80-210-13-04	4.0 mm	2 buttons	International 3-pin
80-210-18-04	4.0 mm	2 buttons	Erbe VIO, ICC, ACC
80-217-05-04	4.0 mm	1 button	KLS Martin
80-220-02-04	4.0 mm	without button, via footswitch	KLS Martin
80-220-41-04	4.0 mm	without button, via footswitch	Bovie
80-220-31-04	4.0 mm	without button, via footswitch	Erbe VIO, ICC, ACC
80-210-21-04	2.4 mm	2 buttons	International 3-pin
80-210-23-04	2.4 mm	Rocker switch	International 3-pin
Electrode handpieces with detachable cord			
80-227-00-04	4.0 mm	2 buttons	-
Connection cord			
80-227-01-04	-	-	KLS Martin
80-227-02-04	-	-	Erbe VIO, ICC, ACC
80-227-03-04	-	-	International 3-pin
Dental handpieces with permanent cord			
80-214-13-04	1.6 mm	1 button	KLS Martin
80-214-12-04	1.6 mm	without button, via footswitch	KLS Martin

The features at a glance



Unbeatable ergonomics, easy handling. See for yourself!

Cut and Coag activated
with pushbuttons
or rocker switches.
Long-life premium quality.

Premium quality materials
- High dielectric strength
- Large number of sterilization cycles
- Highly resistant against cleaning agents

**Anti-kink
grommet**
at cable
transition*

**Premium quality
cord, length 4 m**

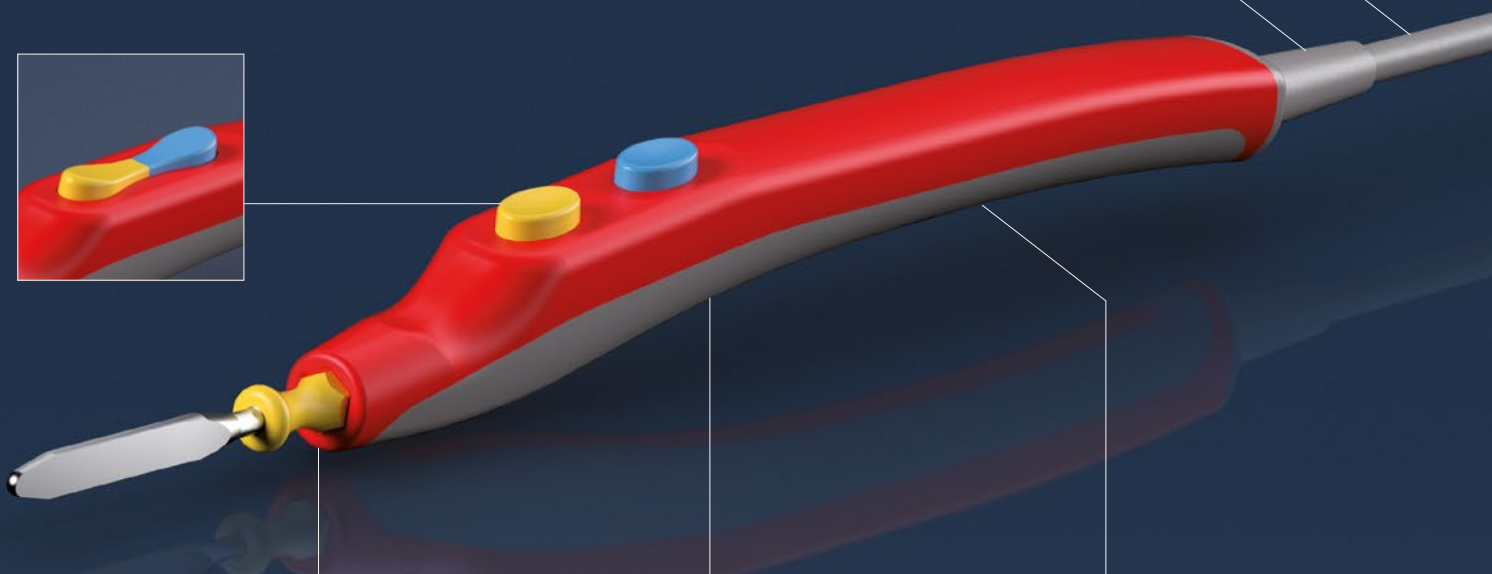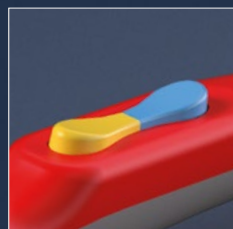


Hex-locking
electrode holder
for electrodes with
Ø 2.4 mm and 4.0 mm

Round electrode holder for
electrodes with Ø 1.6 mm

Premium quality pliant component
So the handpiece will never slip out
of your hand.

**Ergonomically tailored
handpiece**

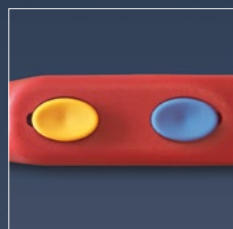

The openings are part of the cleaning concept
and ensure that the smallest amounts of liquid
can be drained off during the cleaning process.

Different plugs for plug
sockets by various device
manufacturers



* in corresponding models the electrode handpiece
may be detached from the cord.

KLS Martin Group

KLS Martin Australia Pty Ltd.

Sydney · Australia
Tel. +61 2 9439 5316
australia@klsmartin.com

KLS Martin do Brasil Ltda.

São Paulo · Brazil
Tel. +55 11 3554 2299
brazil@klsmartin.com

KLS Martin Medical (Shanghai) International Trading Co. Ltd.

Shanghai · China
Tel. +86 21 5820 6251
china@klsmartin.com

KLS Martin India Pvt Ltd.

Chennai · India
Tel. +91 44 66 442 300
india@klsmartin.com

Martin Italia S.r.l.

Milan · Italy
Tel. +39 039 605 67 31
italia@klsmartin.com

Nippon Martin K.K.

Tokyo · Japan
Tel. +81 3 3814 1431
nippon@klsmartin.com

KLS Martin SE Asia Sdn. Bhd.

Penang · Malaysia
Tel. +604 505 7838
malaysia@klsmartin.com

KLS Martin de México S.A. de C.V.

Mexico City · Mexico
mexico@klsmartin.com

Martin Nederland/Marned B.V.

Huizen · The Netherlands
Tel. +31 35 523 45 38
nederland@klsmartin.com

Gebrüder Martin GmbH & Co. KG

Moscow · Russia
Tel. +7 499 792-76-19
russia@klsmartin.com

KLS Martin Taiwan Ltd.

Taipei 106 · Taiwan
Tel. +886 2 2325 3169
taiwan@klsmartin.com

Gebrüder Martin GmbH & Co. KG

Dubai · United Arab Emirates
Tel. +971 4 454 16 55
middleeast@klsmartin.com

KLS Martin UK Ltd.

London · United Kingdom
Tel. +44 1189 000 570
uk@klsmartin.com

KLS Martin LP

Jacksonville · Florida, USA
Tel. +1 904 641 77 46
usa@klsmartin.com

Do you know how to get all important information
about the KLS Martin Energy Devices?

Please download KLS Martin App Energy Devices!
The App is available for Android and iOS.



[https://itunes.apple.com/de/app/
klsmartin-energy-devices/id1198171415?l=de&ls=1&mt=8](https://itunes.apple.com/de/app/klsmartin-energy-devices/id1198171415?l=de&ls=1&mt=8)



[https://play.google.com/store/apps/
details?id=com.klsmartin.energydevices](https://play.google.com/store/apps/details?id=com.klsmartin.energydevices)

Gebrüder Martin GmbH & Co. KG

A company of the KLS Martin Group

KLS Martin Platz 1 · 78532 Tuttlingen · Germany
P.O. Box 60 · 78501 Tuttlingen · Germany
Tel. +49 7461 706-0 · Fax +49 7461 706-193
info@klsmartin.com · www.klsmartin.com

